

# Human VSNL1/VILIP-1 ELISA Kit

Catalog NO.: RK02511

version: 2.0

This package insert must be read in its entirety  
before using this product

## **Background**

VSNL1/ VILIP-1 is a member of the visinin/recoverin subfamily of neuronal calcium sensor proteins. The sequence of which remain same homologous among human , mouse and rat according to GenBank data. VILIP-1 has widespread distributions and high expression levels in many regions containing basal ganglia , procerebrum , cerebellum and cerebrum, and in most neuronal cells of brain.

The encoded protein is strongly expressed in granule cells of the cerebellum where it associates with membranes in a calcium-dependent manner and modulates intracellular signaling pathways of the central nervous system by directly or indirectly regulating the activity of adenylyl cyclase. Alternatively spliced transcripts encoding different proteins have been described. Although VSNL1/ VILIP-1 is constitutively expressed in brain tissues and a few other organs, the expression abundance of VSNL1/ VILIP-1 in brain is significantly higher than that in other parts.

VSNL1/ VILIP-1 plays a complex role in pre-apoptosis and regulation of various signal cascades. VSNL1/ VILIP-1 was constitutively expressed in the brain and a few other organs, but the abundance of VSNL1/ VILIP-1 in the brain was

significantly higher than that in other organs. Therefore, various types of brain diseases, such as ischemic or hemorrhagic, infarction or proliferation, malignant or benign tumors, spontaneous or traumatic, may show abnormal expression of VSNL1/ VILIP-1 (up or down regulation, or phase change).

## **Introduction**

The kit is a sandwich enzyme immunoassay for in vitro quantitative measurement of VSNL1/VILIP-1 in human serum, plasma, cell culture supernates and other biological fluids.

<b>Features</b>	<b>Specifications</b>
Sensitivity	8.625 pg/mL
Detection Range	31.25-2000 pg/mL
SPECIFICITY	There was no obvious cross-reaction or interference between VSNL1/VILIP-1 and its analogues.

## **Principle of the Assay**

This assay employs the quantitative sandwich enzyme immunoassay technique. An antibody specific for human VSNL1/VILIP-1 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any VSNL1/VILIP-1 present is bound by the immobilized antibody. After washing away any unbound substances, and then a detection antibody specific for VSNL1/VILIP-1 is added to the wells and binds to the combination of capture antibody VSNL1/VILIP-1 in sample. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Following incubation and wash steps, a substrate solution is added to the wells and color develops in proportion to the amount of VSNL1/VILIP-1 bound in the initial step. The color development is stopped and the absorbance is measured.

## **Material Provided & Storage Conditions**

Store unopened kit at 2-8 °C for 1 year. Do not use after expiration date. It is highly recommended to use the remaining reagents within 1 month after opening.

<b>Part</b>	<b>Size</b>	<b>Cat. No.</b>	<b>Storage of opened/reconstituted material</b>
Human VSNL1/VILIP-1 Microwell Plate Coated	8x12	RM95136	Put the unused slats back in the aluminum foil bag with the desiccant and reseal them. They can be stored at 2-8°C for 1 month.
Human VSNL1/VILIP-1 Standard Lyophilized	2 vials	RM95137	It is not recommended to use again after redissolving.
Human VSNL1/VILIP-1 Concentrated Biotin Conjugate Antibody (100x)	1 x120µL	RM95138	Store at 2-8°C for 1 month *

Streptavidin-HRP Concentrated (100x)	1 x120μL	RM95139	Store at 2-8°C for 1 month *
Standard/Sample Diluent (R1)	1 x20mL	RM00023	
Biotin-Conjugate Antibody Diluent (R2)	1 x12mL	RM00024	Store at 2-8°C for 1 month *
Streptavidin-HRP Diluent(R3)	1 x12mL	RM00025	
Wash Buffer(20x)	1 x30mL	RM00026	
TMB Substrate	1 x12mL	RM00027	
Stop Solution	1 x6mL	RM00028	
Plate Sealers	4 Strips		
Specification	1		

## **Other Supplies Required**

1. Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 630 nm or 570 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squirt bottle, manifold dispenser, or automated microplate washer.
5. Incubator.
6. Test tubes for dilution of standards and samples.



## **Precautions**

1. Any variation in diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
2. Variations in sample collection, processing, and storage may cause sample value differences.
3. Reagents may be harmful, if ingested, rinse it with an excess amount of tap water.
4. Stop Solution contains strong acid. Wear eye, hand, and face protection.
5. Please perform simple centrifugation to collect the liquid before use.
6. Do not mix or substitute reagents with those from other lots or other sources.
7. Adequate mixing is particularly important for good result. Use a mini-vortexer at the lowest frequency.
8. Mix the sample and all components in the kits adequately, and use clean plastic container to prepare all diluents.
9. Both the sample and standard should be assayed in duplicate, and reagents should be added in sequence in accordance with the requirement of the specification.
10. Reuse of dissolved standard is not recommended.

11. The kit should not be used beyond the expiration date on the kit label.
12. The kit should be away from light when it is stored or incubated.
13. To reduce the likelihood of blood-borne transmission of infectious agents, handle all serum, plasma, and other biological fluids in accordance with NCCLS regulations.
14. To avoid cross contamination, please use disposable pipette tips.
15. Please prepare all the kit components according to the Specification. If the kits will be used several times, please seal the rest strips and preserve with desiccants. Do use up within 2 months.
16. This assay is designed to eliminate interference by other factors present in biological samples.
17. Until all factors have been tested in this assay, the possibility of interference cannot be excluded.
18. The 48T kit is also suitable for the specification.

## **Sample Collection & Storage**

The sample collection and storage conditions listed below are intended as general guidelines. Sample stability has not been evaluated.

Samples containing the correlated IgG as in this kit may interfere with this assay.

**Cell Culture Supernates:** Remove particulates by centrifugation. Assay immediately or aliquot and store samples at  $\leq -20$  °C. Avoid repeated freeze-thaw cycles.

**Serum:** Use a serum separator tube (SST) and allow samples to clot for 30 minutes at room temperature before centrifugation for 15 minutes at 1000xg. Remove serum and assay immediately or aliquot and store samples at  $\leq -20$  °C. Avoid repeated freeze-thaw cycles.

**Plasma:** Collect plasma using EDTA or Heparin as an anticoagulant. Centrifuge for 15 minutes at 1000 × g within 30 minutes of collection. Assay immediately or aliquot and store samples at  $\leq -20$ °C. Avoid repeated freeze-thaw cycles. (Note: Citrate plasma has not been validated for use in this assay.)

**Other biological fluids :** Centrifuge samples for 20 minutes at  $1,000 \times g$ . Collect the supernatants and assay immediately or store samples in aliquot at  $-20^{\circ}\text{C}$  or  $-80^{\circ}\text{C}$  for later use. Avoid repeated freeze-thaw cycles.

**Note:** It is suggested that all samples in one experiment be collected at the same time of the day. Avoid hemolytic and hyperlipidemia sample for serum and plasma.


## **Reagent Preparation**

Bring all reagents to room temperature before use. If crystals have formed in the concentrate, Bring the reagent to room temperature and mix gently until the crystals have completely dissolved.

**Standard** - Reconstitute the Standard Lyophilized with 1.0 mL Standard/Sample Diluent(R1). This reconstitution produces a stock solution of 3500pg/mL. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making

dilutions.

Use the 3500pg/mL standard stock to produce a dilution series (below) with Standard/Sample Diluent(R1). Mix each tube thoroughly and change pipette tips between each transfer (recommended concentration for standard curve: 2000, 1000, 500, 250, 125, 62.5, 31.25, 0pg/mL). Use diluted standards within 60 minutes of preparation.

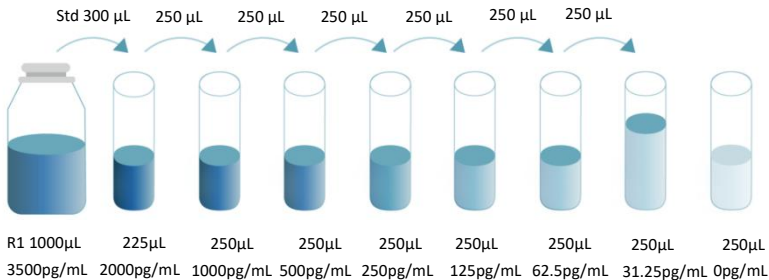

**Working Biotin Conjugate Antibody** - Dilute 1:100 of Concentrated Biotin Conjugate Antibody (100x) with Biotin-Conjugate Antibody Diluent (R2) before use, for example: Add 20 µL of Concentrated Biotin Conjugate Antibody (100x) to 1980 µL Biotin-Conjugate Antibody Diluent (R2) to prepare 2000 µL Working Biotin Conjugate Antibody Buffer.

**Working Streptavidin-HRP** - Dilute 1:100 of Concentrated Streptavidin-HRP (100x) with Streptavidin-HRP Diluent (R3) before use, for example: Add 20  $\mu$  L of Concentrated Streptavidin-HRP (100x) to 1980 $\mu$ L Streptavidin-HRP Diluent (R3) to prepare 2000 $\mu$ L Working Streptavidin-HRP Buffer.

**Wash Buffer** - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Dilute 1:20 with double distilled or deionized water before use, for example: Add 20 mL of Wash Buffer Concentrate to 380 mL of deionized or distilled water to prepare 400 mL of Wash Buffer.

## **Assay Procedure**

Bring all reagents and samples to room temperature before use. It is recommended that all standards, controls, and samples be assayed in duplicate.

1. Prepare all reagents, working standards, and samples as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal properly.
3. Add wash buffer 350  $\mu\text{L}$ /well, aspirate each well after holding 40 seconds, repeating the process two times for a total of three washes.
4. Add 100  $\mu\text{L}$  Standard/sample Diluent (R1) in a blank well.
5. Add 100  $\mu\text{L}$  different concentration of standard or sample in other wells, Cover with the adhesive strip provided. Incubate for 2 hours at 37°C. record the plate layout of standards and sample assay.
6. Prepare the Concentrated Biotin Conjugate Antibody (100x) Working Solution 15 minutes early before use.
7. Repeat the aspiration/wash as in step 3.
8. Add 100  $\mu\text{L}$  Working Biotin Conjugate Antibody in each well, cover with new adhesive Sealer provided. Incubate

for 1 hour at 37°C.

9. Prepare the Streptavidin-HRP Concentrated (100x) Working Solution 15 minutes early before use.
10. Repeat the aspiration/wash as in step 3.
11. Add 100  $\mu$ L Working Streptavidin-HRP in each well, cover with new adhesive Sealer provided. Incubate for 0.5 hour at 37°C.
12. Repeat the aspiration/wash as in step 3.
13. During the incubation, turn on the microplate reader to warm up for 30 minutes before measuring.
14. Add 100  $\mu$ L TMB Substrate to each well. Incubate for 15-20 minutes at 37°C. Protect from light.
15. Add 50  $\mu$ L Stop Solution, determine the optical density of each well within 5 minutes, using a Microplate reader set to 450 nm. If wavelength correction is available, set to 570 nm or 630 nm. If wavelength correction is not available, subtract readings at 570 nm or 630 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may cause higher value and less accurate result.



## **Assay Procedure Summary**

Prepare the standard and reagents

Wash 3 times



Add 100ul of standards or test samples to each well

Incubate for 2 hours at 37°C, then wash 3 times



Add 100ul Working Biotin Conjugate Antibody

Incubate for 1 hour at 37°C, then wash 3 times



Add 100ul Working Streptavidin-HRP

Incubate for 0.5 hour at 37°C, then wash 3 times



Add 100ul Substrate Solution

Incubate for 15-20 min at 37°C under dark condition



Add 50ul Stop Solution




Detect the optical density within 5 minutes under 450nm.

Correction Wavelength set at 570nm or 630nm

## **Calculation of Results**

1. Average the duplicate readings for each standard, control and sample, and subtract the average zero standard optical density (O.D.).
2. Create a standard curve by reducing the data using computer software capable of generating a four-parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the Y-axis against the concentration on the X-axis and draw a best fit curve through the points on a log/log graph. The data may be linearized by plotting the log of the VSNL1/VILIP-1 concentrations versus the log of the O.D. on a linear scale, and the best fit line can be determined by regression analysis.
3. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

## Typical Data

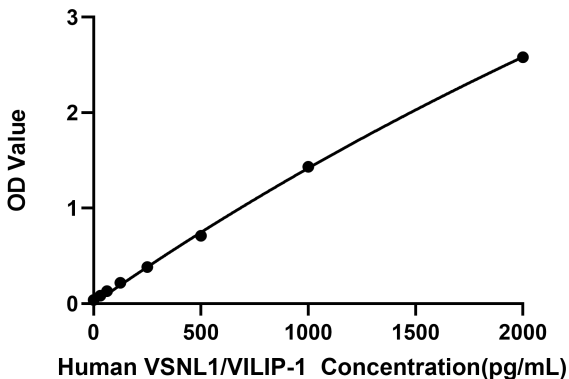

The standard curves are provided for demonstration only. A standard curve should be generated for each set of VSNL1/VILIP-1 assayed.

## **Sensitivity**

The minimum detectable dose (MDD) of VSNL1/VILIP-1 typically less than 8.625pg/mL. The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

## **Specificity**

This method has high sensitivity and specificity for the detection of VSNL1/VILIP-1.

There was no obvious cross-reaction or interference between VSNL1/VILIP-1 and its analogues.

### **Note:**

Due to the limitations of current technology and knowledge, it is impossible to complete the cross-reaction detection between VSNL1/VILIP-1 and all analogues, so the cross-reaction may still exist.

## Precision

### Intra-plate Precision

3 samples with low, middle and high level VSNL1/VILIP-1 were tested 16 times on one plate, respectively.

Intra-Assay: CV<10%

### Inter-plate Precision

3 samples with low, middle and high level VSNL1/VILIP-1 were tested on 3 different plates, 16 replicates in each plate.

Inter-Assay: CV<15%

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	16	16	16	16	16	16
Mean(pg/mL)	92	353	1391	75	325	1406
Standard deviation	0.01	0.04	0.14	0.14	0.41	1.22
CV(%)	6.1	7.8	7.6	7.5	8.8	7.0

## **Recovery**

Matrices listed below were spiked with certain level of VSNL1/VILIP-1 and the recovery rates were calculated by comparing the measured value to the expected amount of VSNL1/VILIP-1 in samples.

Sample	Average Recovery (%)	Range (%)
Cell Culture Media(n=5)	93	88-106
Serum(n=5)	87	81-102

## Linearity

Five samples were spiked with known concentrations of target proteins, and the recovery experiments were performed to obtain the range of recovery rates and average recovery rates. Five samples were diluted 2, 4, 8 and 16 times for recovery experiments, and the range of recovery rates and average recovery rates were obtained.

/	/	Cell Culture Media(n=5)	Serum(n=5)
1:2	Average of Expected (%)	89	89
	Range (%)	87-90	82-95
1:4	Average of Expected (%)	102	93
	Range (%)	93-111	91-94
1:8	Average of Expected (%)	99	106
	Range (%)	96-102	105-107
1:16	Average of Expected (%)	96	102
	Range (%)	90-102	98-106



## Trouble Shooting

<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
High Background	Insufficient washing	Sufficiently wash plates as required. Ensure appropriate duration and number of washes. Ensure appropriate volume of wash buffer in each well.
	Incorrect incubation procedure	Check whether the duration and temperature of incubation are set up as required.
	Cross-contamination of samples and reagents	Be careful of the operations that could cause cross-contamination. Use fresh reagents and repeat the tests.
No signal or weak signal	Incorrect use of reagents	Check the concentration and dilution ratio of reagents. Make sure to use reagents in proper order.
	Incorrect use of microplate reader	Warm the reader up before use. Make sure to set up appropriate main wavelength and correction wavelength.
	Insufficient colour reaction time	Optimum duration of colour reaction should be limited to 15-25 minutes.

	Read too late after stopping the colour reaction	Read the plate in 5 minutes after stopping the reaction.
	Matrix effect of samples	Use positive control.
Too much signal	Contamination of TMB substrate	Check if TMB substrate solution turns blue. Use new TMB substrate solution.
	Plate sealers reused	Use a fresh new sealer in each step of experiments.
	Protein concentration in sample is too high	Do pre-test and dilute samples in optimum dilution ratio.
Poor Duplicates	Uneven addition of samples	Check the pipette. Periodically calibrate the pipette.
	Impurities and precipitates in samples	Centrifuge samples before use.
	Inadequate mixing of reagents	Mix all samples and reagents well before loading.

\*For research purposes only. Not for therapeutic or diagnostic purposes.